


***Instruction Manual and
Experiment Guide for
the PASCO scientific
Model ME-8569***

012-04774A
1/92

DENSITY SET

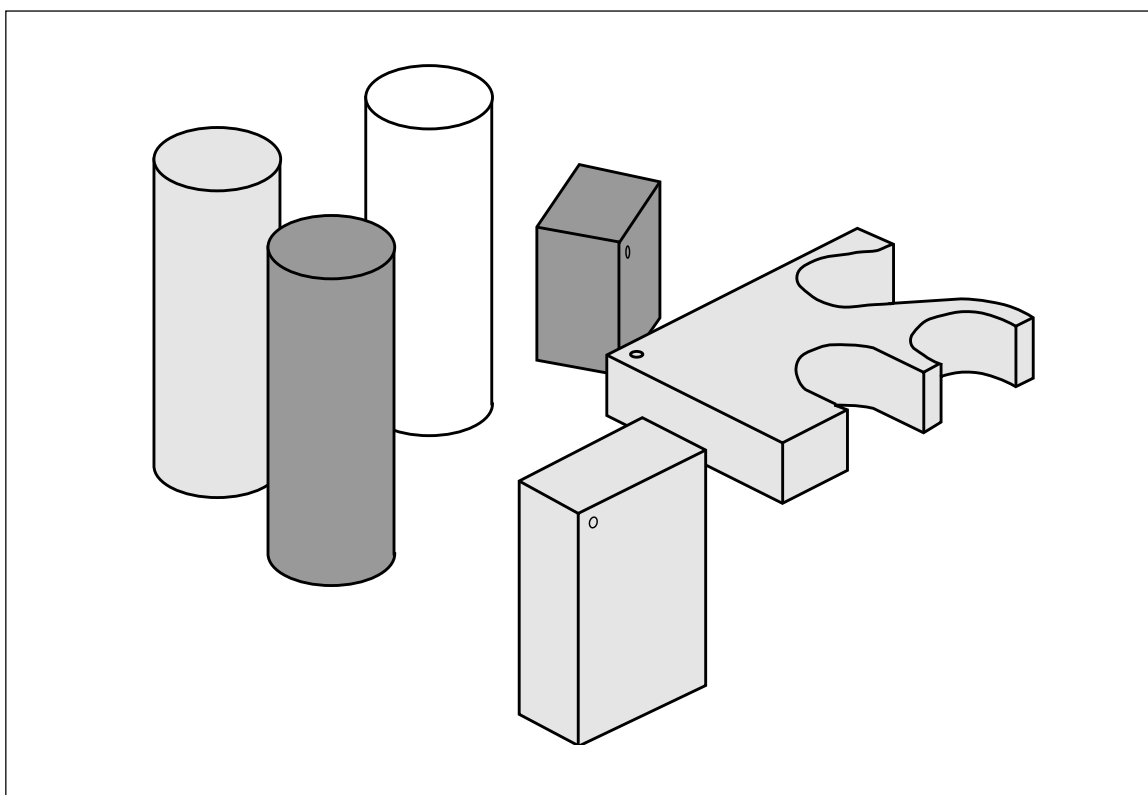

Copyright © January 1992

\$5.00

Table of Contents

Section	Page
Copyright, Warranty, and Equipment Return	ii
Introduction	1
Equipment	1
Replacement Parts	2
Experiments	
Experiment 1: Length, Area, and Volume	3
Experiment 2: Density	5
Experiment 3: Buoyancy Force	7
Experiment 4: Specific Heat	11

Copyright and Warranty

Please—Feel free to duplicate this manual subject to the copyright restrictions below.

Copyright Notice

The PASCO scientific Model ME-8569 Density Set manual is copyrighted and all rights reserved. However, permission is granted to non-profit educational institutions for reproduction of any part of this manual providing the reproductions are used only for their laboratories and are not sold for profit. Reproduction under any other circumstances, without the written consent of PASCO scientific, is prohibited.

Limited Warranty

PASCO scientific warrants this product to be free from defects in materials and workmanship for a period of one year from the date of shipment to the customer. PASCO

will repair or replace, at its option, any part of the product which is deemed to be defective in material or workmanship. This warranty does not cover damage to the product caused by abuse or improper use. Determination of whether a product failure is the result of a manufacturing defect or improper use by the customer shall be made solely by PASCO scientific. Responsibility for the return of equipment for warranty repair belongs to the customer. Equipment must be properly packed to prevent damage and shipped postage or freight prepaid. (Damage caused by improper packing of the equipment for return shipment will not be covered by the warranty.) Shipping costs for returning the equipment, after repair, will be paid by PASCO scientific.

Equipment Return

Should this product have to be returned to PASCO scientific, for whatever reason, notify PASCO scientific by letter or phone BEFORE returning the product. Upon notification, the return authorization and shipping instructions will be promptly issued.

NOTE: NO EQUIPMENT WILL BE ACCEPTED FOR RETURN WITHOUT AN AUTHORIZATION.

When returning equipment for repair, the units must be packed properly. Carriers will not accept responsibility for damage caused by improper packing. To be certain the unit will not be damaged in shipment, observe the following rules:

1. The carton must be strong enough for the item shipped.
2. Make certain there is at least two inches of packing material between any point on the apparatus and the inside walls of the carton.
3. Make certain that the packing material can not shift in the box, or become compressed, thus letting the instrument come in contact with the edge of the box.

Introduction

The PASCO Model ME-8569 Density Set allows useful experiments in density, buoyancy force, and specific heat as well as providing an effective means of distinguishing between length, volume, and area.

Equipment

The density set consists of six objects:

- aluminum cylinder
- aluminum block
- brass cylinder
- brass block
- polypropylene cylinder
- aluminum irregular shape

Each object has a small hole through which a string can be tied to suspend the object. The objects in this set can be categorized in three different ways: those objects having the same volume, same mass, or same density.

SAME VOLUME

- aluminum cylinder
- aluminum block
- brass cylinder
- polypropylene cylinder
- aluminum irregular shape

SAME MASS

- aluminum cylinder
- aluminum block
- brass block
- aluminum irregular shape

SAME DENSITY

- aluminum cylinder
- aluminum block
- aluminum irregular shape

SAME DENSITY

- brass cylinder
- brass block

Replacement Parts

<u>ITEM</u>	<u>PASCO PART NUMBER</u>
aluminum cylinder	648-04768
aluminum block	648-04772
brass cylinder	648-04770
brass block	648-04771
polypropylene cylinder	648-04769
aluminum irregular shape	648-04773

Experiment #1: Length, Area, and Volume

REQUIRED EQUIPMENT:

Density set
 Calipers (SF-8711)
 String (SE-8050)
 Overflow can (SE-8568)
 Beaker for catching water
 Graduated cylinder (50 ml)

PURPOSE:

The purpose of this experiment is to distinguish between length, area, and volume.

PROCEDURE

I. LENGTH

Using the calipers, measure the longest side of each of the five regularly-shaped objects. Record the results in Table 1.1.

II. AREA

A. Using the calipers, measure the diameter of each of the three cylinders. Divide the diameter by two to get the radius, r .

Calculate the area of the circular end of the cylinders using:

$$A = \pi r^2.$$

Record the results in Table 1.1.

B. Using the calipers, measure the width and height of one end of each of the blocks. Calculate the area of the end of the block by multiplying the width by the height. Record the results in Table 1.1.

III. VOLUME

A. By Calculation

For each of the regularly-shaped objects, calculate the volume by multiplying the area of one end by the length of the object. Record in Table 1.1.

B. By Displacement of Water

For each of the regularly-shaped objects, find the volume by finding the volume of water that each one displaces:

1. Put the beaker under the overflow can spout as shown in Figure 1.
2. Pour water into the overflow can until it overflows into the beaker. Allow the water to stop overflowing on its own and empty the beaker into the sink and return it to its position under the overflow can spout without jarring the overflow can.
3. Tie a string to each of the objects (including the irregularly-shaped object).
4. Gently lower the first object into the overflow can until it is completely submerged. Allow the water to stop overflowing and then pour the water from the beaker into the graduated cylinder.

Table 1.1: Length, Area, and Volume

OBJECT	LENGTH	AREA	VOLUME	DISPLACED VOLUME	% DIFF.
AL. CYL.					
AL. BLOCK					
BRASS CYL.					
BRASS BLOCK					
PLASTIC					
AL. IRREG. SHAPE					

Measure the volume of water that was displaced by reading the water level in the graduated cylinder in milliliters ($1 \text{ ml} = 1 \text{ cm}^3$). Record this volume in Table 1.1.

5. Repeat this procedure for the other objects. Note that the plastic object will float in water so it cannot be used in this part of the experiment.

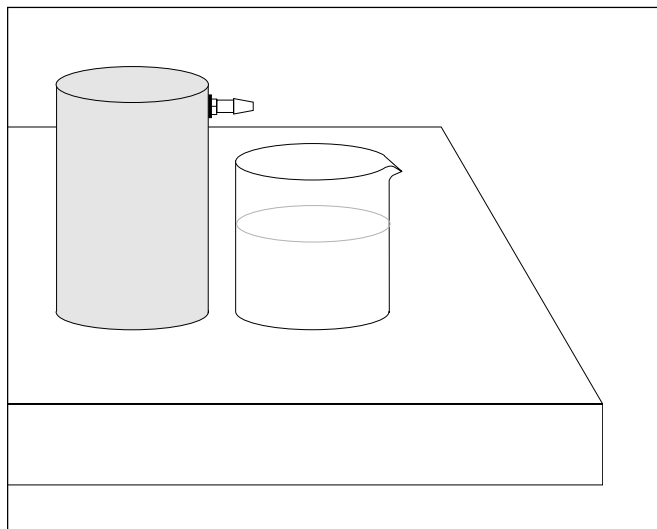


Figure 1: Overflow Setup

ANALYSIS

For each of the regular objects, calculate the percent difference between the two values found for the volume.

QUESTIONS

1. Which objects have nearly the same volume?
2. Which objects have nearly the same length?
3. Which objects have nearly the same cross-sectional area?
4. Did any two objects have the same volume but did not have the same length or the same area?

Experiment #2: Density

REQUIRED EQUIPMENT:

Density set: brass cylinder and brass block

Calipers (SF-8711)

String (SE-8050)

Overflow can (SE-8568)

Beaker for catching water

Graduated cylinder (50 ml)

Triple-beam balance (SE-8708)

PURPOSE:

This experiment finds the densities of different-shaped objects made of the same material.

PROCEDURE

Using the triple-beam balance, find the mass of the brass cylinder and the brass block. Record the results in Table 2.1.

I. CALCULATED VOLUME

A. Cylinder

Using the calipers, measure the length, L , of the brass cylinder. Record the results in Table 2.1. Measure the diameter of the cylinder. Divide the diameter by two to get the radius, r . Calculate the volume of the cylinder using:

$$V = \pi r^2 L.$$

Record the results in Table 2.1.

B. Block

Using the calipers, measure the length (L), width (W), and height (H) of the brass block. Calculate the volume of the block using:

$$V = L \times W \times H.$$

Record the results in Table 2.1.

II. DISPLACED VOLUME

For each of the two brass objects, find the volume by finding the volume of water that each one displaces:

1. Put the beaker under the overflow can spout as shown in Figure 1.
2. Pour water into the overflow can until it overflows into the beaker. Allow the water to stop overflowing on its own and empty the beaker into the sink and return it to its position under the overflow can spout without jarring the overflow can.
3. Tie a string to each of the objects.
4. Gently lower the first object into the overflow can until it is completely submerged. Allow the water to stop overflowing and then pour the water from the beaker into the graduated cylinder.

Measure the volume of water that was displaced by reading the water level in the graduated cylinder in milliliters ($1 \text{ ml} = 1 \text{ cm}^3$). Record this volume in Table 2.1.

5. Repeat this procedure for the other object.

ANALYSIS

For each method of finding the volume, calculate the density, d , of each object using:

$$d = \frac{\text{mass}}{\text{volume}}$$

Obtain the accepted value for the density of brass from your instructor or from a reference book.

For each experimental value of the density, calculate the percent difference from the accepted value and record the results in Table 2.2.

Table 2.1: Measurements

OBJECT	MASS	DIMENSIONS	VOLUME	DISPLACED VOLUME
CYLINDER				
BLOCK				

QUESTIONS

1. Which object took up more space?
2. Which object weighed more?
3. Which object was more dense?

Table 2.2: Results

ACCEPTED DENSITY VALUE = _____

OBJECT	CALCULATED DENSITY	DENSITY BY DISPLACEMENT	% DIFFERENCE FROM ACCEPTED
CYLINDER			
BLOCK			

Experiment #3: Bouyant Force

REQUIRED EQUIPMENT:

Density set

String (SE-8050)

Overflow can (SE-8568)

Beaker for catching water

Graduated cylinder (50 ml)

Triple-beam balance (SE-8708)

PURPOSE:

This experiment measures the buoyant force on an object in water by using Archimedes' Principle and by finding the upward force on the object while it is submerged.

PROCEDURE

I. ARCHIMEDES' PRINCIPLE

Archimedes' Principle states that the buoyant force exerted on an object partially or fully submerged in a fluid will be equal to the weight of the fluid displaced by the object. To use this principle to find the buoyant force exerted on each object, follow these steps:

1. Find the mass of the empty beaker and record this at the top of Table 3.1. Put the beaker under the overflow can spout as shown in Figure 1.

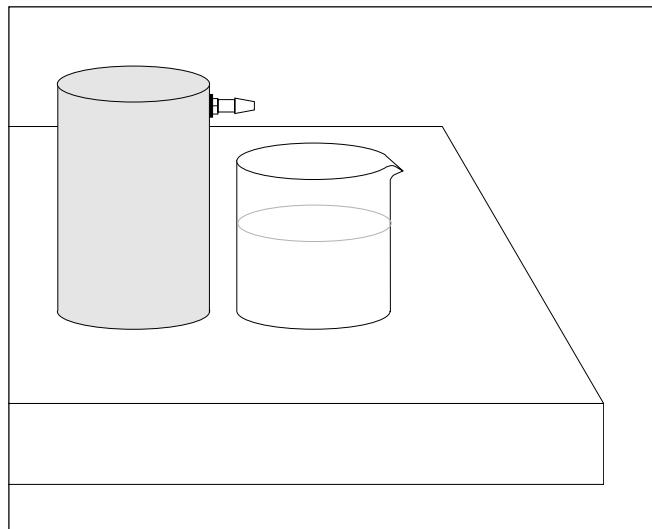


Figure 1: Overflow Setup

2. Pour water into the overflow can until it overflows into the beaker. Allow the water to stop overflowing on its own and empty the beaker into the sink and return it to its position under the overflow can spout without jarring the overflow can.
3. Tie a string to each of the objects.
4. Gently lower the first object into the overflow can until it is completely submerged. Allow the water to stop overflowing and then find the mass of the beaker with the water in it. Record this mass in Table 3.1.

Table 3.1: Measurements

MASS OF BEAKER = _____

OBJECT	MASS OF BEAKER + WATER	MASS OF WATER	WEIGHT OF DISPLACED WATER
AL. CYLINDER			
AL. BLOCK			
BRASS CYL.			
BRASS BLOCK			
PLASTIC CYL.			
AL. IRREG. SHAPE			
BRASS CYL. 1/2 SUBMERGED			

5. Calculate the mass of the water by subtracting the mass of the beaker. Then calculate the weight of the displaced water by multiplying by the acceleration due to gravity (9.8 m/s^2). Record the result in Table 3.1 and Table 3.2.

6. Repeat this procedure for the other objects. Note that the plastic object will float in water but it can still be used in this

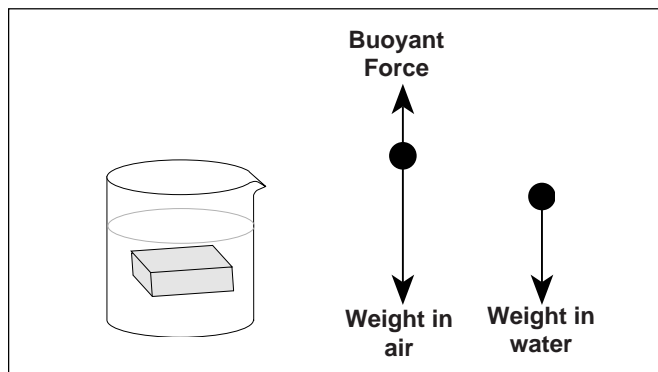


Figure 2: Free Body Diagram

part of the experiment. Also, repeat the procedure for the brass cylinder with only half the cylinder submerged.

II. UPWARD FORCE

When an object is submerged in a fluid, the apparent weight of the object is less than the weight in air because of the buoyant force (See Figure 2).

Thus the buoyant force can be calculated by finding the difference between the weight of the object in air and the apparent weight of the object when it is submerged in water.

1. Put the triple-beam balance on top of a stand as shown in Figure 3. Tie a string to the bottom of the pan.
2. Hang the first object from the string. The balance will read the same as when the object is placed on top of the pan. Multiply the mass by the acceleration due to gravity (9.8 m/s^2) and record the object's weight in Table 3.2.

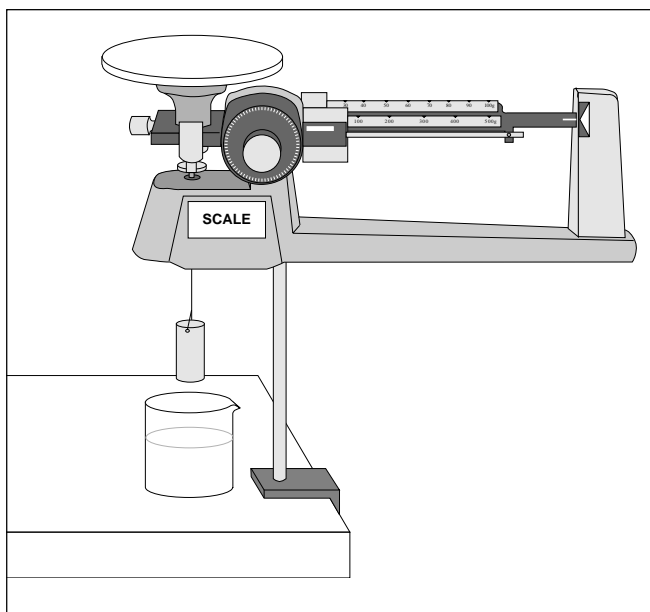


Figure 3: Scale Setup

Table 3.2: Measurements

OBJECT	WEIGHT IN AIR	WEIGHT IN WATER	BF = $W_{\text{AIR}} - W_{\text{WATER}}$	WEIGHT OF DISPLACED WATER
AL. CYLINDER				
AL. BLOCK				
BRASS CYL.				
BRASS BLOCK				
PLASTIC CYL.				
AL. IRREG. SHAPE				
BRASS CYL. 1/2 SUBMERGED				

3. While the object is still hanging from the balance, submerge the object in a beaker of water so that the entire object is under water but it is not touching the sides or bottom of the beaker. Record the reading on the scale and multiply by gravity to get the apparent weight. Record in Table 3.2.
4. Calculate the buoyant force by taking the difference between the weight in air and the weight in water. Record in Table 3.2.
5. Repeat these steps for all the objects. Note that the plastic cylinder will float so don't try to completely submerge it in the water. Also, for the half-submerged brass cylinder, find the apparent weight in the water when only half the cylinder is submerged.

QUESTIONS

1. In each case, is the buoyant force found using the difference between weights equal to the weight of the water displaced?
2. Which objects had the same buoyant force when submerged?
Why?
3. For the plastic cylinder, what was the weight in water?
4. How was the buoyant force for the totally submerged brass cylinder related to the buoyant force for the half-submerged brass cylinder?

Notes

Experiment #4: Specific Heat

EQUIPMENT REQUIRED:

Density set: Aluminum and Brass cylinders and blocks

String (SE-8050)

Beaker for heating water

Styrofoam cup

Cold water

Triple-beam balance (SE-8708)

Bunsen burner

Heating stand

Thermometer (SE-9083)

PURPOSE:

This experiment shows that the specific heat of a material depends on the type of material but not on the amount of material.

PROCEDURE

The specific heat, c , of a material is defined to be the amount of heat needed to raise the temperature of one gram of the material one degree Celsius. To measure this heat, the method of mixtures will be used. The following procedure will be repeated for four objects (aluminum cylinder, aluminum block, brass cylinder, and brass block):

1. Heat a beaker of water to boiling. While the water continues to boil, hang the object by a string so that the object is completely submerged in the boiling water but it is not touching the bottom of the beaker. Allow the object to come to equilibrium with the boiling water (wait about 5 minutes). See Figure 1.

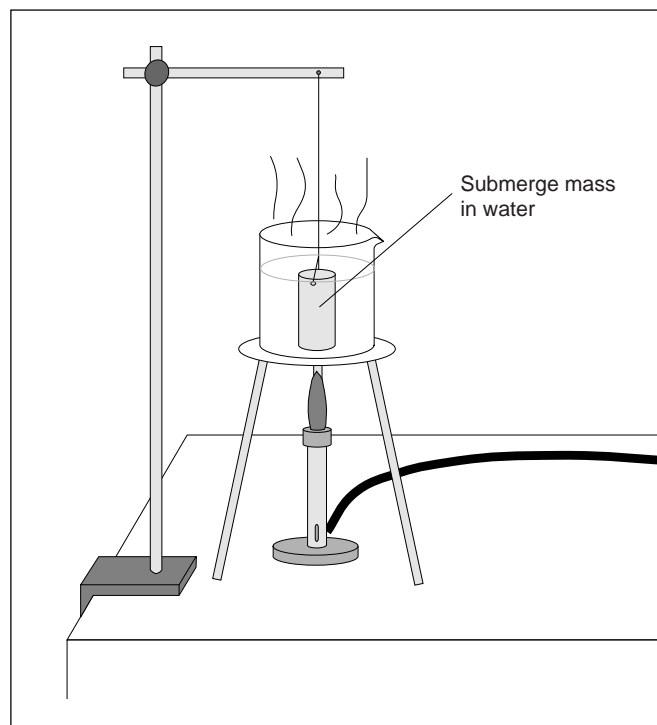


Figure 1: Heating Setup

2. While waiting, find the mass of a dry styrofoam cup and then prepare a styrofoam cup of cold water (about 3°C below room temperature) into which the heated object can be completely submerged. When the object has finished heating, record the temperature of the cold water to the nearest tenth of a degree, then record the temperature of the boiling water, and quickly transfer the heated object from the hot water to the cold water (See Figure 2). Record in Table 4.1.

Table 4.1: Data

MASS OF CUP = _____

OBJECT	HOT TEMPERATURE	COLD TEMPERATURE	EQUILIBRIUM TEMPERATURE	MASS OF WATER
AL. CYL.				
AL. BLOCK				
BRASS CYL.				
BRASS BLOCK				

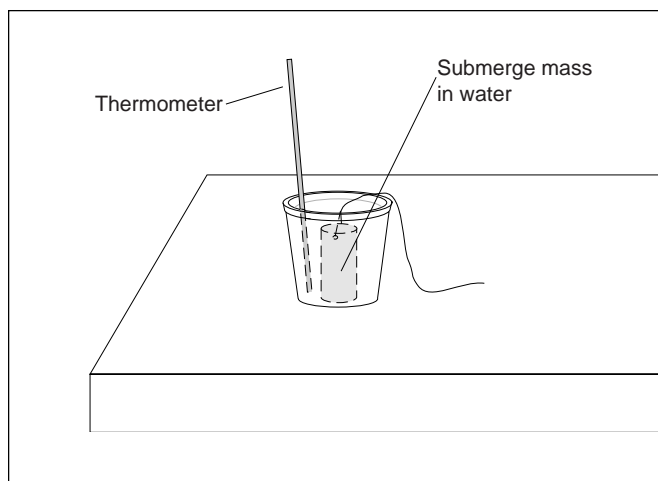


Figure 2: Equilibrium Temperature

3. Watch the temperature of the cold water rise, stirring the water gently. After several minutes, the temperature will peak out (above room temperature) and begin to fall. Record this equilibrium (peak) temperature to the nearest tenth of a degree in Table 1.

4. After the equilibrium temperature has been reached, remove the object from the water and weigh the water with the cup. Subtract the mass of the cup from the mass of the water with the cup and record the mass of the water in Table 4.1.

5. Calculate the change in temperature for the object:

$$\Delta T = \text{HOT TEMP} - \text{EQUILIBRIUM TEMP.}$$

Record this in Table 4.2.

6. Calculate the change in temperature for the water in the styrofoam cup:

$$\Delta T = \text{EQUILIBRIUM TEMP} - \text{COLD TEMP.}$$

Record this in Table 4.2.

7. When the object is cooled in the water, the heat lost by the object is equal to the heat gained by the water because energy is conserved (neglecting any losses to the environment):

$$\Delta Q_{\text{OBJECT}} = \Delta Q_{\text{WATER}}$$

$$(mc\Delta T)_{\text{OBJECT}} = (mc\Delta T)_{\text{WATER}}$$

Solving for the specific heat of the object gives:

$$c = \frac{m_{\text{WATER}} c_{\text{WATER}} (\Delta T)_{\text{WATER}}}{m_{\text{OBJECT}} (\Delta T)_{\text{OBJECT}}}$$

where c_{WATER} is 1 cal/g °C. Calculate the specific heat for the object and record the results in Table 4.2.

8. Look up the accepted value for the specific heat for the material used and record in Table 4.2. Calculate the percent difference between the experimental value and the accepted value.

9. Repeat this procedure for the other objects.

QUESTIONS

1. Do the two aluminum objects have the same specific heat?
2. Do the two brass objects have the same specific heat?
3. How does starting the cold water below room temperature minimize the effect of the losses to the environment?

Table 4.2: Calculations

OBJECT	OBJECT ΔT	WATER ΔT	SPECIFIC HEAT	ACCEPTED VALUE	% DIFF.
AL. CYL.					
AL. BLOCK					
BRASS CYL.					
BRASS BLOCK					

Technical Support

Feedback

If you have any comments about the product or manual, please let us know. If you have any suggestions on alternate experiments or find a problem in the manual, please tell us. PASCO appreciates any customer feedback. Your input helps us evaluate and improve our product.

To Reach PASCO

For technical support, call us at 1-800-772-8700 (toll-free within the U.S.) or (916) 786-3800.

fax: (916) 786-3292

e-mail: techsupp@PASCO.com

web: www.pasco.com

Contacting Technical Support

Before you call the PASCO Technical Support staff, it would be helpful to prepare the following information:

➤ If your problem is with the PASCO apparatus, note:

- Title and model number (usually listed on the label);
- Approximate age of apparatus;
- A detailed description of the problem/sequence of events. (In case you can't call PASCO right away, you won't lose valuable data.);
- If possible, have the apparatus within reach when calling to facilitate description of individual parts.

➤ If your problem relates to the instruction manual, note:

- Part number and revision (listed by month and year on the front cover);
- Have the manual at hand to discuss your questions.

